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This cooling bag keeps the liquid in your medication between 2 and 8 degrees with PCM cooling elements from the freezer. No power, no battery and no maintenance on the go. Read this manual in full once before your first trip; after that, packing takes just a few minutes.

## Your model at a glance

| Model           | Gel packs | Cooling time | Capacity                          |
|-----------------|-----------|--------------|-----------------------------------|
| LifeinaBag12    | 1         | 12 hours     | 3 injection pens                  |
| iCool Weekender | 1         | 12 hours     | 3 injection pens                  |
| LifeinaBag24    | 2         | 24 hours     | 5 injection pens                  |
| iCool Prestige  | 2         | 24 hours     | 5 injection pens                  |
| iCool VacPack   | 3         | 18 hours     | inner dimensions 180 x 70 x 70 mm |

All models maintain 2–8°C with the bag closed for the stated cooling time, measured on liquid temperature at 24°C outside temperature.

## Ready to go in four steps

### 1 Freeze the gel packs

Place all gel packs flat in the freezer at minus 18 degrees for at least 12 hours. After 12 hours they are fully charged. Longer is fine: store them there permanently, then they are always ready for use.

### 2 Check and condition

Not sure your freezer is at minus 18 degrees, or is there frost or ice on the elements? Place the gel packs next to the cooling bag in the open air for 5 to 30 minutes until the frost has gone. Do this too after more than 12 hours in the freezer: longer freezing sometimes forms ice on the surface. This way the bag starts at the right temperature and also protects your medication against freezing.

### 3 Pack the bag

Place the gel packs in the designated compartment and take your medication straight from the refrigerator (around 5 degrees). Per model:

- LifeinaBag12 and iCool Weekender: the gel pack in the zippered compartment.
- LifeinaBag24 and iCool Prestige: two gel packs. One in the silver-coloured compartment, one in the flap on the inside, where the zipper is.
- iCool VacPack: one gel pack on the bottom, two against the sides and inside them the blue tray with your medication.

### 4 Close the bag and keep it closed

Zip the bag fully closed and keep it closed as much as possible on the way. The insulation does the work as long as the bag stays shut; every opening lets warm air in and shortens the cooling time. Do not overfill the bag: tension on the zipper pushes it open a fraction and cold air escapes through that gap.

## | How to keep your bag cool the longest

- Keep the bag out of direct sun and away from heat sources such as the parcel shelf, the heating or the car bonnet.
- The higher the outside temperature, the harder the gel packs work. On a hot day you only get the maximum cooling time with a bag that stays closed and in the shade.
- Open the bag only when you need your medication and close it again right after.
- Simply keep the cooling bag in your backpack on the way. The backpack buffers the ambient temperature and helps keep the temperature in the bag low for longer.
- At your destination? Put your medication in the fridge and the gel packs in the freezer. The next morning you repack and set off fully charged again.

Good to know: the stated cooling time was tested at an outside temperature of 24 degrees, the average temperature in an aircraft cabin. Every degree of difference means about 30 minutes of cooling time. At 30 degrees outside the difference is 6 degrees, so the bag keeps your medication cool for about 3 hours less. If the outside temperature is lower, the bag cools even longer.

## | PCM cooling elements: cooling without freezing risk

Every cooling bag in this manual comes with PCM cooling elements as standard. PCM stands for phase change material: a material with a freezing point of 3 degrees. The element therefore freezes and thaws around 3 degrees and releases its cold exactly in the range where your medication belongs. As a result the bag stays stable between 2 and 8 degrees for a long time and your medication runs no risk of freezing. The instructions for freezing and use are in the step-by-step plan on page 2.

Bag and elements belong together. Every cooling bag is precisely matched to the supplied cooling elements and the stated amount of medication. Using extra elements, other elements or our elements in another bag? That combination has not been tested and we make no claims about the cooling time or temperature range.

## | Safety in case of damage or leakage

The gel in the PCM elements has very low toxicity and is harmless to your health in normal use. If an element does get damaged, follow these instructions.


- Skin contact: wash the skin immediately with water and soap.
- Eye contact: rinse thoroughly with plenty of water.
- Swallowed: rinse the mouth well with water and drink plenty of water. Consult a doctor if discomfort persists.

### Waste disposal and the environment

- Never rinse the gel down the sink or the toilet. The gel contains thickening agents that swell and can seriously block pipes or the sewer.
- Wipe up spilled gel with paper or a cloth. Careful: the gel makes the surface slippery right away.
- Dispose of the broken element and the cleaned-up remains with the general waste.
- Always replace a damaged or leaking element before your next trip. Order separate PCM elements at [www.medicooler.nl](http://www.medicooler.nl).

## Measuring temperature: air and liquid

The cooling bag is designed to keep the liquid in your medication between 2 and 8 degrees, not the air in the bag. Air reacts instantly to a frozen gel pack or an opened zipper and therefore fluctuates considerably. Liquid reacts much more slowly due to its greater thermal mass. In the graph you see the difference over the full cooling time: the liquid drops to 2 degrees in the first hour, then stays stable for a long time and only rises faster towards 8 degrees near the end. The air even shoots below freezing right after packing, then climbs much faster and ends well above 8 degrees, while your medication stays neatly within range the whole time.

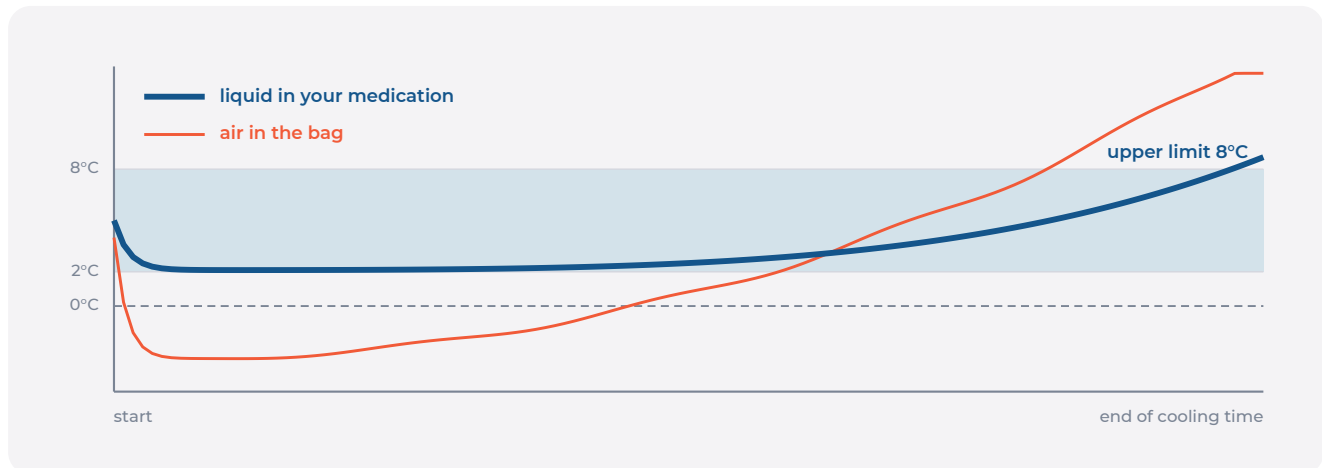If you measure with a loose air sensor in the bag, you will see the orange line. That says little about your medication itself. Want to test the bag? Place the sensor in a small bottle filled with water, close the bag and leave it closed for 24 hours. That way you measure what counts: the temperature of the liquid. Official test reports and comparison graphs are available via our customer service.

## Travel, flying and security screening

- All MediCooler cooling bags are allowed as hand luggage. Always keep your medication with you in the cabin, never in the hold.
- Frozen gel packs are not liquids at security screening and are always allowed in your hand luggage.
- Your medication and cooling elements are medically necessary and therefore fall under the exemption to the 100 ml rule. Show them separately at screening if asked.
- For the medication itself customs sometimes asks for a doctor's letter; for the cooling bag it is never needed. With your order you will find the MediCooler medical certificate for your doctor to fill in and stamp.
- Travelling to a destination outside Europe? Check the country's rules for bringing medication in advance.
- Travelling back or onward without using the cooling bag? Then pack the gel packs in your checked luggage. When not frozen they count as liquids and are not allowed in hand luggage.

## Condensation on the bag

- Condensation can form on the outside of the bag. The more air flows past the bag, the more condensation: warm air that cools down against the cold bag condenses into moisture.
- Want to catch the moisture? Put the cooling bag in a plastic bag.

## | Multi-day trips

The routine for every travel day: put your medication in the fridge of your accommodation in the evening and the gel packs in the freezer. The next morning the elements are recharged and you repack the bag in a few minutes. No freezer at your accommodation? Ask reception; almost every hotel freezes gel packs for guests. This way you extend the cooling time of your bag without limit, day after day.

## | Maintenance and storage

- Clean the bag only with a damp cloth. Never wash the bag by hand or in the washing machine.
- Never iron the bag and never dry it on the heating or near another heat source. Just let it air dry.
- Let the bag dry fully open after use and after cleaning before storing it. This prevents damp spots and odours.
- Store the bag open or loosely closed, in a dry place out of the sun.
- Keep the gel packs in the freezer, then they are always ready for departure. Outside the freezer, store them at room temperature, away from sharp objects.
- Check the elements for small cracks or leaks before every trip and replace them when in doubt.

## | Good to know: what affects the cooling time

We designed and tested these cooling bags with the greatest care. Even so, there are circumstances that affect the performance of every passive cooling bag and over which we as a manufacturer have no control. Keep in mind:

- the outside temperature and humidity
- the number of times the bag is opened
- how long the bag stays open
- the temperature of the gel packs when packing
- the temperature of the medication when packing
- exposure to direct sunlight
- how long the medication stays in the bag
- the position of gel packs and medication

### DISCLAIMER

The responsibility for keeping medication at the correct temperature lies entirely with the user. MediCooler is in no case liable for direct or indirect damage to property, medication or health arising from the use or incorrect use of the cooling bag or the cooling elements. Doubts about your purchase after receiving it? Then contact our customer service at [info@medicooler.nl](mailto:info@medicooler.nl) before first use to return the bag.

## | Service and warranty

All MediCooler products are covered by the statutory warranty. Questions about your trip, your medication or your cooling bag? Email [info@medicooler.nl](mailto:info@medicooler.nl). You will get a personal answer from a specialist.

We wish you a good trip. Your medication stays at the right temperature.

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